

SCANLAN® Sundt™

Rotatable & Non-Rotatable
Graduated Suction System




INTERNATIONAL

Surgical products designed and manufactured by the Scanlan family since 1921

SCANLAN® Sundt™ Rotatable & Non-Rotatable Graduated Suction System

The SCANLAN® SUNDT™ Graduated Suction System provides complete control to atraumatically retract, probe, lift, and evacuate fluid or tissue from the surgical field. The SUNDT™ system offers a choice of seven tip diameters, each available in three lengths. The universal suction handle easily accommodates all tip diameters, lengths, and configurations.


Universal Suction Handle (shown actual size)

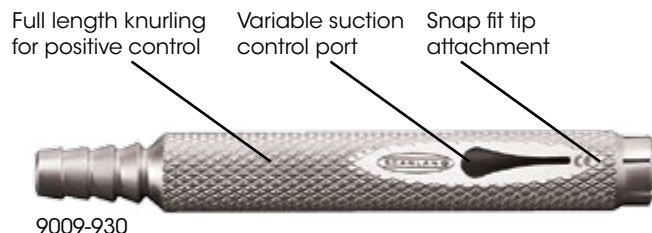

Features

- Atraumatic tips may be used to retract or probe delicate tissue
- Precise fingertip control for atraumatic tissue manipulation or aspiration
- Each suction tip is tapered to reduce clogging
- One universal suction handle accommodates every tip configuration
- Tips and handle easily separate for cleaning and sterilization
- Seven suction tip diameters, ranging from 3 Fr. to 11 Fr.
- Three suction tip lengths
- Rotating or Non-rotating tips
- Guide pin securely holds tips in place at angle (for non-rotating tips)
- Universal SCANLAN® V.I.P.™ Sterilization tray holds 7 tips plus handles, fits all lengths
- Tips, handle, and tray may be purchased separately or in sets

Three Suction Tip Lengths (shown actual size)

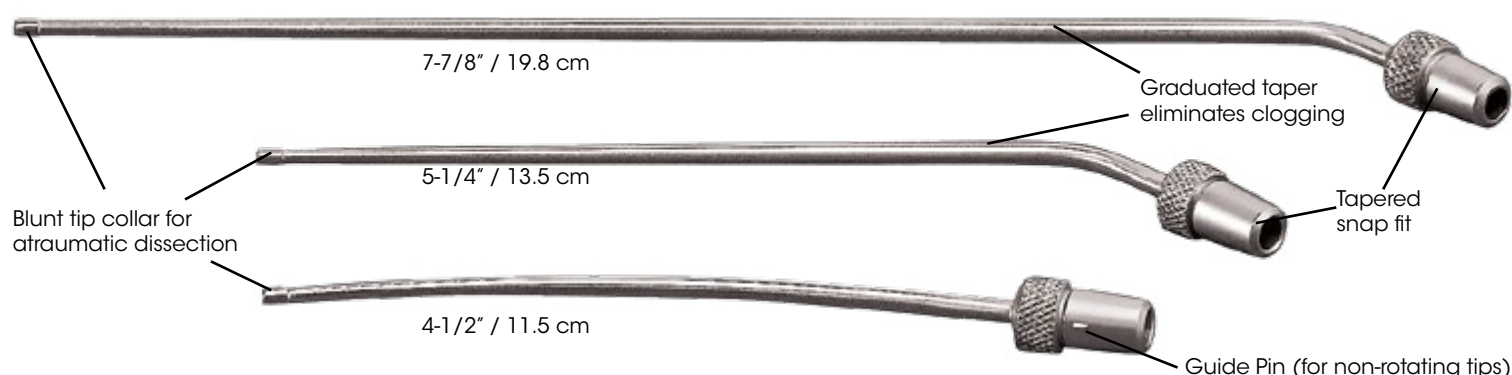

The SCANLAN® Sundt™ Graduated Suction System may be purchased as a set or as individual components.	Rotatable Sets			Non-Rotatable Sets		
	Each set has all seven suction tips and includes one universal suction handle and one V.I.P.™ sterilization / storage case.			Each set has all seven suction tips and includes one universal suction handle and one V.I.P.™ sterilization / storage case.		
	9009-940	9009-960	9009-980	9009-950	9009-970	9009-990
SCANLAN® Sundt™ Components	4-1/2" / 11.5 cm	5-1/4" / 13.5 cm	7-7/8" / 19.8 cm	4-1/2" / 11.5 cm	5-1/4" / 13.5 cm	7-7/8" / 19.8 cm
• 3 Fr. Tip	9009-941	9009-961	9009-981	9009-951	9009-971	9009-991
• 4 Fr. Tip	9009-942	9009-962	9009-982	9009-952	9009-972	9009-992
• 5 Fr. Tip	9009-943	9009-963	9009-983	9009-953	9009-973	9009-993
• 6 Fr. Tip	9009-944	9009-964	9009-984	9009-954	9009-974	9009-994
• 7 Fr. Tip	9009-945	9009-965	9009-985	9009-955	9009-975	9009-995
• 9 Fr. Tip	9009-947	9009-967	9009-987	9009-957	9009-977	9009-997
• 11 Fr. Tip	9009-949	9009-969	9009-989	9009-959	9009-979	9009-999
Universal Suction Handle (3" / 7.5 cm) fits all tips	9009-930	9009-930	9009-930	9009-930	9009-930	9009-930
SCANLAN® Sundt™ V.I.P.™ Sterilization / Storage Case 6" x 10" x 1.5" / 15.2 x 25.4 x 3.8 cm, holds 7 tips plus handles, fits all lengths	9009-934	9009-934	9009-934	9009-934	9009-934	9009-934



INTERNATIONAL

Surgical products designed and manufactured by the Scanlan family since 1921
 ONE SCANLAN PLAZA • SAINT PAUL, MINNESOTA 55107 • USA
 800-328-9458 • 651-298-0997 • FAX 651-298-0018 • www.scanlaninternational.com

© 2016 Scanlan International, Inc.
 MK-B-N-SUNDT REV. C